






Universitat de Barcelona




Unitat de RMN

## Descripción del modelo

### Tema-1b





Universitat de Barcelona

## Remembering

An atom is classified according to its number of protons and neutrons:

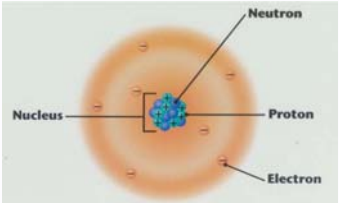
**The number of protons** in the nucleus of an atom determines an element's atomic number "chemical element"

All carbon atoms, and only carbon atoms, contain six protons and have an atomic number of 6

**The number of neutrons** determines the "isotope" of that element. All atoms have a mass number, which is the sum of protons and neutrons.

The carbon atom has several isotopes. The most abundant with six neutrons and one with seven neutrons ( $^{12}\text{C}$  and  $^{13}\text{C}$ )

The nuclei of all atoms may be characterized by:  
a nuclear spin quantum number ( $I$ )



Only nuclei with ***spin number ( $I$ )  $\neq 0$***  can absorb/emit electromagnetic radiation.

Unitat de RMN Centres Científics i Tecnològics

## NMR Active nucleus

Universitat de Barcelona

---

The nuclear spin quantum number  $I$  can either be equal to zero, or to multiples of  $1/2$

For atoms with there is no nuclear spin and therefore, they cannot have a nuclear magnetic resonance. These atoms are called NMR silent. All other values of  $I$  yield nuclear spin.

**Even mass # and Even atomic #**  
 $I=0$  ( $^{12}\text{C}$ ,  $^{16}\text{O}$ , etc )

**Odd mass# and Even atomic #**  
 $I=1/2$  ( $^1\text{H}$ ,  $^{13}\text{C}$ ,  $^{15}\text{N}$ )

**Odd mass# and Odd atomic#**  
 $I= n/2$  ( $n \neq 1$ ) integer  $I=3/2$  ( $^{11}\text{B}$ ,  $^{23}\text{Na}$ );  $^{53}\text{Ca}$   $I=7/2$

**Even mass # and Odd atomic #**  
 $I=\text{whole integer}$   $I=1$  ( $^2\text{H}$ ,  $^{14}\text{N}$ ) ;  $I=3$   $^{10}\text{B}$

No Nuclear spin  
**NMR Inactive**

Nuclear Spin  
 Spherical charge distribution

Ellipsoidal charge distribution

Unitat de RMN Centres Científics i Tecnològics

## Tabla periódica

Universitat de Barcelona

---

**Tabla periódica de isótopos de RMN.**

Número de masa  $\rightarrow$  ← Número atómico  
 ← Símbolo químico

1																	18
$^1\text{H}$																	$^4\text{He}$
$^7\text{Li}$	$^9\text{Be}$															$^{19}\text{F}$	$^{20}\text{Ne}$
$^{23}\text{Na}$	$^{24}\text{Mg}$	$^{27}\text{Al}$	$^{28}\text{Si}$	$^{29}\text{P}$	$^{32}\text{S}$	$^{35}\text{Cl}$	$^{36}\text{Ar}$									$^{79}\text{Br}$	$^{80}\text{Kr}$
$^{39}\text{K}$	$^{40}\text{Ca}$	$^{45}\text{Sc}$	$^{48}\text{Ti}$	$^{51}\text{V}$	$^{52}\text{Cr}$	$^{55}\text{Mn}$	$^{56}\text{Fe}$	$^{59}\text{Co}$	$^{63}\text{Ni}$	$^{65}\text{Cu}$	$^{68}\text{Zn}$	$^{70}\text{Ga}$	$^{72}\text{Ge}$	$^{75}\text{As}$	$^{78}\text{Se}$	$^{81}\text{Br}$	$^{84}\text{Kr}$
$^{85}\text{Rb}$	$^{87}\text{Sr}$	$^{88}\text{Y}$	$^{90}\text{Zr}$	$^{91}\text{Nb}$	$^{92}\text{Mo}$	$^{93}\text{Tc}$	$^{98}\text{Ru}$	$^{101}\text{Rh}$	$^{102}\text{Pd}$	$^{106}\text{Ag}$	$^{108}\text{Cd}$	$^{112}\text{In}$	$^{114}\text{Sn}$	$^{115}\text{Sb}$	$^{117}\text{Te}$	$^{126}\text{I}$	$^{129}\text{Xe}$
$^{133}\text{Cs}$	$^{137}\text{Ba}$	$^{138}\text{La}$	$^{140}\text{Ce}$	$^{141}\text{Pr}$	$^{146}\text{Nd}$	$^{150}\text{Pm}$	$^{151}\text{Sm}$	$^{157}\text{Eu}$	$^{162}\text{Gd}$	$^{167}\text{Tb}$	$^{172}\text{Dy}$	$^{175}\text{Ho}$	$^{177}\text{Er}$	$^{178}\text{Tm}$	$^{173}\text{Yb}$	$^{175}\text{Lu}$	
$^{223}\text{Fr}$	$^{227}\text{Ra}$	$^{227}\text{Ac}$	$^{232}\text{Th}$	$^{232}\text{Pa}$	$^{235}\text{U}$	$^{237}\text{Np}$	$^{241}\text{Pu}$	$^{243}\text{Am}$	$^{247}\text{Cm}$	$^{251}\text{Bk}$	$^{255}\text{Cf}$	$^{259}\text{Es}$	$^{263}\text{Fm}$	$^{267}\text{Md}$	$^{271}\text{No}$	$^{275}\text{Lr}$	


Espin nuclear =  $I$

1/2	1	3/2	5/2	3	7/2	4	9/2	5	6	7
-----	---	-----	-----	---	-----	---	-----	---	---	---

Armando Ariza Castolo aariza@cinvestav.mx

Unitat de RMN Centres Científics i Tecnològics

## Nuclear properties related to NMR


  
 Universitat de Barcelona


Nuclear properties

Nuclide	Spin $I$	Electric quadrupole moment <sup>a)</sup> [eQ] [10 <sup>-28</sup> m <sup>2</sup> ]	Natural abundance [%]	Relative sensitivity <sup>b)</sup>	Gyromagnetic ratio $\gamma^{a)}$ [10 <sup>7</sup> rad T <sup>-1</sup> s <sup>-1</sup> ]	NMR frequency [MHz] <sup>b)</sup> ( $B_0 = 2.3488$ T)
<sup>1</sup> H	1/2	—	99.985	1.00	26.7519	100.0
<sup>2</sup> H	1	2.87 x 10 <sup>-3</sup>	0.015	9.65 x 10 <sup>3</sup>	4.1066	15.351
<sup>3</sup> H <sup>c)</sup>	1/2	—	—	1.21	28.5350	106.664
<sup>6</sup> Li	1	-6.4 x 10 <sup>-4</sup>	7.42	8.5 x 10 <sup>-3</sup>	3.9371	14.716
<sup>10</sup> B	3	8.5 x 10 <sup>-2</sup>	19.58	1.99 x 10 <sup>-2</sup>	2.8747	10.746
<sup>11</sup> B	3/2	4.1 x 10 <sup>-2</sup>	80.42	0.17	8.5847	32.084
<sup>12</sup> C	0	—	98.9	—	—	—
<sup>13</sup> C	1/2	—	1.108	1.59 x 10 <sup>-2</sup>	6.7283	25.144
<sup>14</sup> N	1	1.67 x 10 <sup>-2</sup>	99.63	1.01 x 10 <sup>-3</sup>	1.9338	7.224
<sup>15</sup> N	1/2	—	0.37	1.04 x 10 <sup>-3</sup>	-2.7126	10.133
<sup>16</sup> O	0	—	99.96	—	—	—
<sup>17</sup> O	5/2	-2.6 x 10 <sup>-2</sup>	0.037	2.91 x 10 <sup>-2</sup>	-3.6280	13.557
<sup>19</sup> F	1/2	—	100	0.83	25.1815	94.077
<sup>23</sup> Na	3/2	0.1	100	9.25 x 10 <sup>-2</sup>	7.0704	26.451
<sup>24</sup> Mg	5/2	0.22	10.13	2.67 x 10 <sup>-3</sup>	-1.6389	6.1195
<sup>28</sup> Si	1/2	—	4.70	7.84 x 10 <sup>-3</sup>	-5.3190	19.865
<sup>31</sup> P	1/2	—	100	6.63 x 10 <sup>-2</sup>	10.8394	40.481
<sup>39</sup> K	3/2	5.5 x 10 <sup>-2</sup>	93.1	5.08 x 10 <sup>-4</sup>	1.2499	4.667
<sup>43</sup> Ca	7/2	-5.0 x 10 <sup>-2</sup>	0.145	6.40 x 10 <sup>-3</sup>	-1.8028	6.728
<sup>57</sup> Fe	1/2	—	2.19	3.37 x 10 <sup>-5</sup>	0.8687	3.231
<sup>59</sup> Co	7/2	0.42	100	0.28	6.3015	23.614
<sup>119</sup> Sn	1/2	—	8.58	5.18 x 10 <sup>-2</sup>	-10.0318	37.272
<sup>133</sup> Cs	7/2	-3.0 x 10 <sup>-3</sup>	100	4.74 x 10 <sup>-2</sup>	3.5339	13.117
<sup>195</sup> Pt	1/2	—	33.8	9.94 x 10 <sup>-3</sup>	5.8383	21.499

**NMR Periodic Table**  
[http://www-usr.rider.edu/~grushow/nmr/NMR\\_tutor/periodic\\_table/nmr\\_pt\\_frameset.html](http://www-usr.rider.edu/~grushow/nmr/NMR_tutor/periodic_table/nmr_pt_frameset.html)

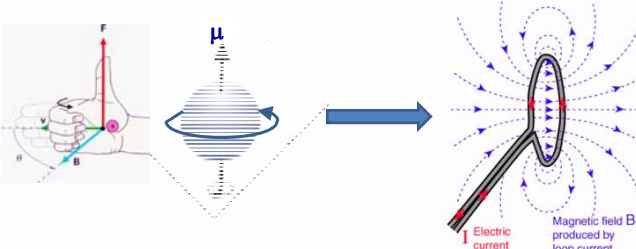
Unitat de RMN Centres Científics i Tecnològics

## Nuclear spin and Magnetic Moment


  
 Universitat de Barcelona

Nucleus rotates about its axis (spin)  
 Nuclear spin results in angular momentum (**P**).  $P = \hbar \times \sqrt{I(I+1)}$   
*Quantized spin quantum number I*  
*2I+1 States I, I-1, I-2, ..., -I*

Since the nucleus is charged, spin will produce a magnetic momentum (**u**)  
 $\mu = \gamma P$  Where  $\gamma$  is the proportionality constant called the gyromagnetic ratio




Similar to magnetic field created by electric current flowing in a coil

↑ Electric current  
 Magnetic field B produced by loop current

The magnetic moment is quantized (m)  
 for proton  $m = +\frac{1}{2}$  &  $-\frac{1}{2}$   $m = I, I-1, I-2, \dots, -I$

Unitat de RMN Centres Científics i Tecnològics

## Quantization of magnetic moment



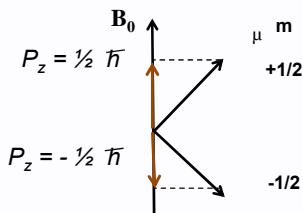
The angular momentum for a nucleus, in a static magnetic field, will be oriented directionally

$$P_z = m\hbar$$

$m$  is the magnetic quantum number  $m = I, I-1, I-2, \dots, -I$

therefore, there are  $(2I+1)$  values for  $m$ , and  $(2I+1)$  possible orientations for the angular momentum

For a nucleus with  $I = 1/2$  there are two possible orientations of the magnetic moment




This is the simplest case only two orientations and only one transition between them

the components of the magnetic moment along z

$$\mu_z = \gamma P_z = m\gamma\hbar$$


Unitat de RMN Centres Científics i Tecnològics

## Precession Classical description



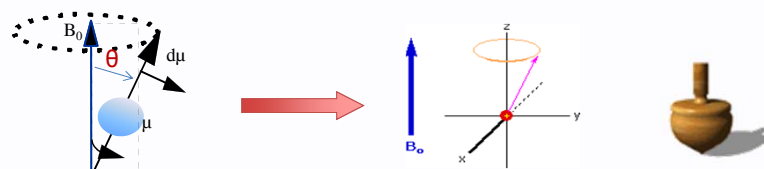In an external magnetic field the nuclear spin precesses at **Larmor frequency**

$$\omega_0 = -\gamma B_0 \text{ (rad s}^{-1}\text{)} \xrightarrow{\nu = \omega_0 / 2\pi} \nu_0 = \gamma B_0 / 2\pi \text{ (Hz)}$$




**Larmor frequency related with  $\gamma$  and  $B_0$**

The relative orientation of the magnetic moment ( $\theta$ ) depends on the value of  $I$

For a nucleus with  $I = 1/2$  there are two possible orientations of the magnetic moment


**(General rule: # orientations =  $2I+1$ ):**

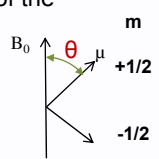This is the simplest case, only two orientations and only one transition between them




Unitat de RMN Centres Científics i Tecnològics


## Magnetic alignment



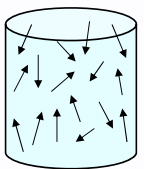
$\hbar = \gamma h / 4\pi$



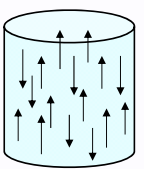In the absence of external field, each nuclei is energetically degenerate



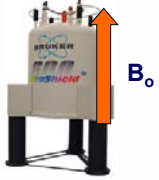Add a strong external field ( $B_0$ ) and the nuclear magnetic moment: aligns with (low energy) against (high-energy)




$B_0$

At equilibrium (in the magnetic field), there is excess of nuclei in the a state of low energy

The energy of a spin in a magnetic field (E) will depend on a static magnetic field called  $B_0$ , and  $\mu$ .

$$E = -\mu_z B_0 = -m\gamma\hbar B_0$$

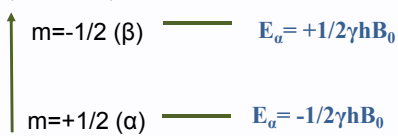For  $l=1/2$  then  $m$  can be  $+1/2$  and  $-1/2$




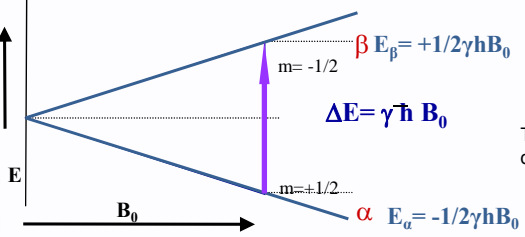
$m = -1/2 (\beta) \quad E_\beta = +1/2\gamma\hbar B_0$

$m = +1/2 (\alpha) \quad E_\alpha = -1/2\gamma\hbar B_0$

Unitat de RMN Centres Científics i Tecnològics

## Nuclear Energy levels in a Static Magnetic Field





$E = -\mu_z B_0 = -m\gamma\hbar B_0$

The energy difference is linearly dependent on  $\gamma$  and on  $B_0$

- At thermal equilibrium the energy difference between  $\alpha$  and  $\beta$  states prevents these states from being equally populated
- The relative population of a particular state is given by the Boltzman distribution:

$$\frac{N_\beta}{N_\alpha} \approx 1 - \frac{\gamma \hbar B_0}{K_B T}$$

$B_0 = 1.41\text{T (60 MHz)} \quad \frac{N_\beta}{N_\alpha} = 0.9999904$

$B_0 = 1.41\text{T (800 MHz)} \quad \frac{N_\beta}{N_\alpha} = 0.999987$

$N_\alpha$  = number of nuclei in the  $\alpha$  state

$N_\beta$  = number of nuclei in the  $\beta$  state

$K_B$  = Boltzman constant  $K_B = 1.3085 \times 10^{-23} \text{ JK}^{-1}$


$T$  = temperature (Kelvin)

The Difference is very small (ppm)

**NMR IS INSENSITIVE**

Unitat de RMN Centres Científics i Tecnològics

## Net Magnetization

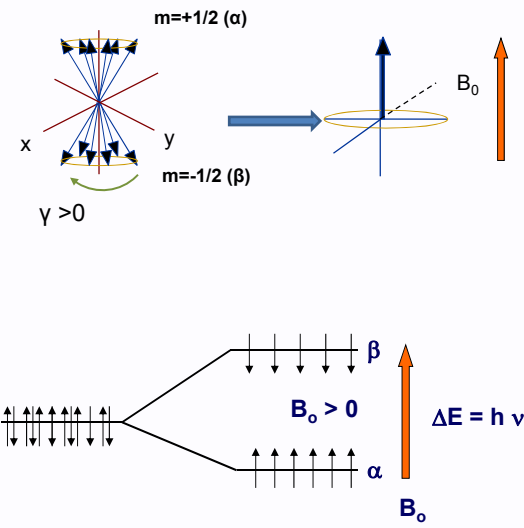


Classic vectorial View:

- Nuclei either align with or against external magnetic field along the z-axis.
- Since more nuclei align with field, net magnetization ( $M_0$ ) exists parallel to external magnetic field


Quantum Description:

- Nuclei either populate low energy ( $\alpha$ , aligned with field) or high energy ( $\beta$ , aligned against field)
- Net population in a lower energy level.
- Absorption of radio-frequency promotes nuclear spins from  $\alpha \rightarrow \beta$ .



Unitat de RMN Centres Científics i Tecnològics

## Resonance



The relationship between the Larmor frequency and  $\Delta E$  is as follows:

$$\Delta E = \hbar \nu_i \quad \Delta E = \gamma \hbar B_0 \quad \text{and} \quad \nu_L = (\gamma/2\pi) B_0, \quad \Delta E = h \nu_L$$

Transitions between energy levels (between  $\alpha$  and  $\beta$  spin states for spin 1/2 nuclei) are quantized, and can only be promoted by an energy  $\Delta E$

In NMR the transitions are promoted by an applied electromagnetic field,  $B_1$  with a frequency,  $\nu_i$  (radiofrequency) matching the Larmor frequency of the nucleus

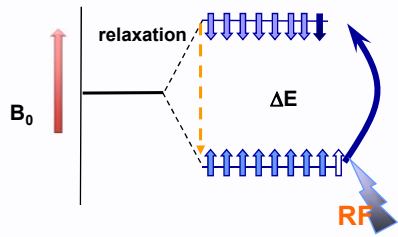
$$\nu_L = \nu_i = (\gamma/2\pi) B_0 \quad \Delta E = h \nu_i$$

This is known as **resonance** ( $\nu_i = \nu_L$ ), when the frequency of our externally applied electromagnetic field ( $B_1$ ) are coincident with the Larmor frequency of the nucleus of interest

$N_\alpha > N_\beta$  →  $\nu_i = \nu_L$  signal

$N_\alpha = N_\beta$  → Saturation

<http://ochem.jsd.claremont.edu/movies.dir/nmr.htm>



Unitat de RMN Centres Científics i Tecnològics




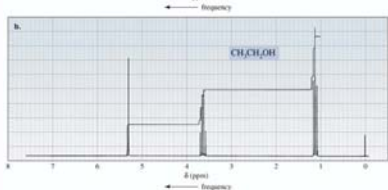
## Continuous Wave (CW)

U  
B Universitat de Barcelona

NMR can be performed like other spectroscopy's (UV/vis, IR) by simply slowly varying the frequency of monochromatic incident radiation and monitoring for absorption

The first spectrometers using a frequency sweep, a constant magnetic field, to obtain the spectra.

One alternative was to sweep the magnetic field while maintaining a constant frequency.

The scan rate determines the resolution and the possibility of saturation of the signals. The average duration of a sweep might be between 2 and 10 minutes.

The system is very inefficient because much of the time was recorded only noise in the spectrum.

Most spectrometers could not add different sweeps

Unitat de RMN Centres Científics i Tecnològics

## Fourier Transform

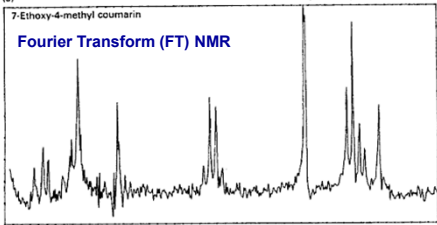
U  
B Universitat de Barcelona

How to efficiently detect a range (spectrum) of NMR frequencies

4.3 SENSITIVITY OF FOURIER SPECTROSCOPY 157

CH<sub>3</sub>-CH<sub>2</sub>O-  
CCOC1=CC=C(C=C1)OC(C)=O (0.01 Molar)  
 Me

(a) 7-Ethoxy-4-methyl coumarin  
 500 Impulse responses of 1 s length  
**Fourier Transform (FT) NMR**



(b) Continuous Wave (CW) NMR  
 1 Scan in 500 s

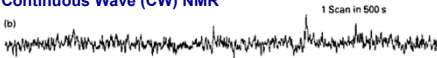



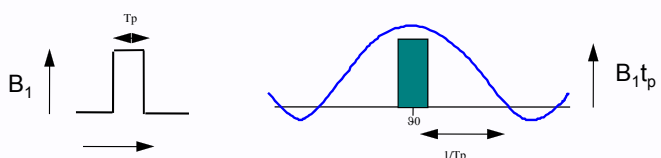
FIG. 4.3.4. 60-MHz proton magnetic resonance spectra of 7-ethoxy-4-methyl coumarin. (a) Fourier transform of 500 free induction signals recorded in 500 s. (b) Single scan recorded in 500 s by slow passage on the same instrument. (Reproduced from Ref. 4.130.)

R.R Ernst et al Principles NMR  
 Unitat de RMN Centres Científics i Tecnològics

## Pulse NMR



In pulsed Fourier transform NMR, for a given nucleus ( $^1\text{H}$  for example), all frequencies are excited simultaneously by a short, high power radiofrequency pulse ( $B_1$  field)



The pulse is applied at a particular frequency,  $\nu_1$ , but a short pulse excites a large continuous band of frequencies (the **bandwidth**) centered around  $\nu_1$ .

The useful or effective bandwidth is proportional to  $1/\tau_p$  ( $\tau_p$  is the **pulse length**, also called the **pulse width or pulse duration**).


The **pulse amplitude** is a measure of the power with which the pulse is applied, and determines the strength of the  $B_1$  field

$\tau_p$  is usually very short, i.e.  $\mu\text{s}$

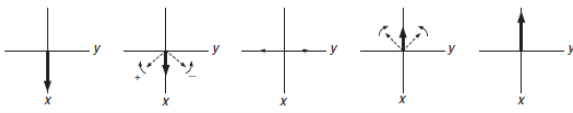
For all frequencies can be excited in a homogeneous mode must be satisfied that  $\gamma B_1 = 2\pi SW$ , (SW equal to spectral window),  $PW_{90} \ll 1/4sw$

Unitat de RMN Centres Científics i Tecnològics

## Rotating Frame




The RF pulse provides an oscillating magnetic field  $B_1$  (at Larmor Frequency) in transverse plane (is equivalent to two counter-rotating vectors)




**Simplification:** Laboratory frame to Rotating Frame

The system rotate to the equal RF frequency  $\nu_0$

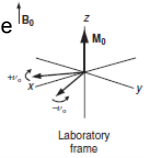The  $-2\nu_0$  is far to the resonance frequency and may be ignored




Laboratory frame

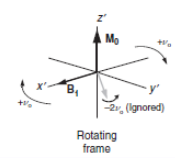


Rotating frame

**Analogy**

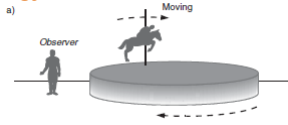

a) Moving

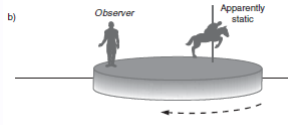Laboratory frame




b) Apparently static

Rotating frame


Unitat de RMN Centres Científics i Tecnològics



## RMN Experiment




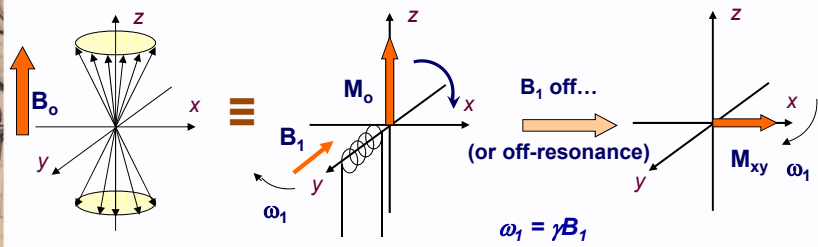**Resonant condition:** frequency ( $\omega_1$ ) of  $B_1$  matches Larmor frequency ( $\omega_0$ ) energy is absorbed and population of  $\alpha$  and  $\beta$  states are perturbed.




$\omega_1 = \gamma B_1$

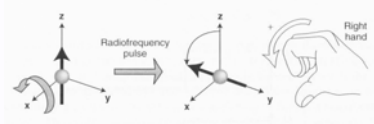$M_0$  now precesses about  $B_1$  for as long as the  $B_1$  field is applied.




Right-hand rule

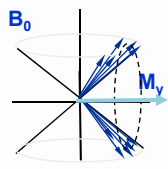

Phase coherence


Unitat de RMN Centres Científics i Tecnològics

## Pulse effect




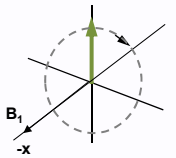The pulse generates  $M_{xy}$  transverse magnetization

The flip angle is:  $\alpha = 360 \gamma B_1 t_p$  degrees



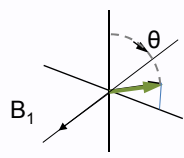
The 90° pulse correspond to maximum signal




The 180° not produce signal (saturation)

Unitat de RMN Centres Científics i Tecnològics

## Pulse phase

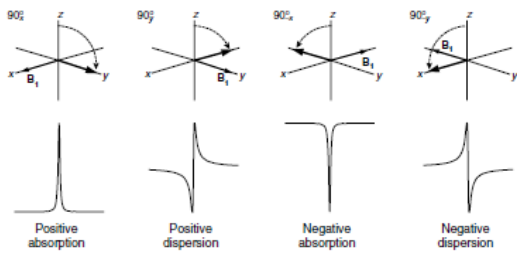
 Universitat de Barcelona

Detection in -y




Positive absorption      Positive dispersion      Negative absorption      Negative dispersion

Some sequences use a different phase pulse: Phase Cycling  
 Selecting the some signals in NMR experiment and rejecting the those that are not required

*Unitat de RMN Centres Científics i Tecnològics*


## NMR Signal Detection
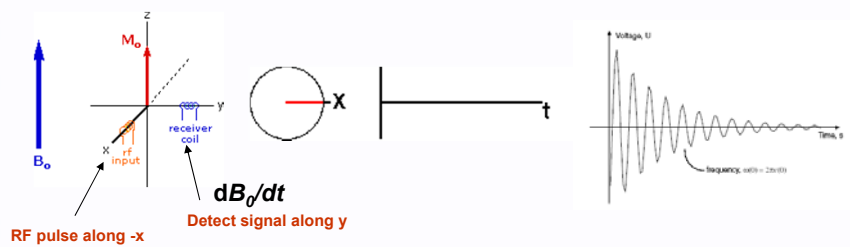
### FID: Free Induction Decay


 Universitat de Barcelona

The pulse generates  $M_{xy}$  transverse magnetization that precesses around the z axis at precession frequency  $\omega_0$

The FID reflects the change in the magnitude of  $M_{xy}$  as the signal is changing relative to the receiver along the y-axis



**RF pulse along -x**      **Detect signal along y**

*Again, the signal is precessing about  $B_0$  at its Larmor Frequency ( $\omega_0$ ).*

*Unitat de RMN Centres Científics i Tecnològics*

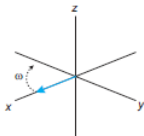
## Events after the pulse: Evolution of the magnetization

The magnetic moment ( $M$ ) precessing about a static magnetic field ( $B_0$ ) results in a local magnetic field ( $B'$ ) varying in time ( $dM/dt \propto dB'/dt$ )


The X & Y components  
(Free Induction Decay):

$$M_y(t) = -M_0 \cos(\omega_0 t) e^{-t/T_2}$$


$$M_x(t) = M_0 \sin(\omega_0 t) e^{-t/T_2}$$



$M_x = M_0 \cos \omega t e^{-t/T_2}$

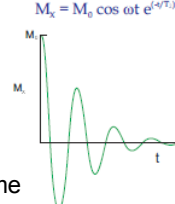

$M_y = M_0 \sin \omega t e^{-t/T_2}$

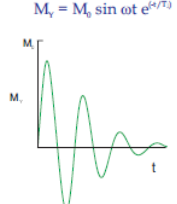

$T_2$  transverse relaxation time

Unitat de RMN Centres Científics i Tecnològics

## NMR Sensitivity in NMR

The sensitivity of a nucleus depends of :

- Gyromagnetic constant
- External magnetic field
- Natural abundance isotope to observe

$$dM/dt \propto \gamma B_0 M \propto N \gamma^3 B_0^2 h^2 I(I+1) / (3k_B T)$$

Isotope	I	$\gamma$ ( $10^7 \text{ rad T}^{-1} \text{ s}^{-1}$ )	Abundancia N(%)	Resonance Frec. B=2.3488T	relative* sensitivity
1H	1/2	26.7519	99.98	100.0	1
19F	1/2	25.1815	100	94.077	0.83
31P	1/2	10.8394	100	40.481	$6.63 \times 10^{-2}$
13C	1/2	6.7283	1.10	25.144	$1.56 \times 10^{-2}$
2H	1	4.1066	0.015	15.351	$9.65 \times 10^{-3}$
15N	1/2	-2.7126	0.37	10.133	$1.04 \times 10^{-3}$

$\gamma \text{ } ^1\text{H} = 26,753 \text{ rad/G}$

$\gamma \text{ } ^{13}\text{C} = 6,728 \text{ rad/G}$

**Ratio  $(\gamma \text{ } ^1\text{H} / \gamma \text{ } ^{13}\text{C})^3 \approx 64$**


If we consider the term A (Natural abundance)      $^1\text{H} \approx 100\%$  ;  $^{13}\text{C} \approx 1\%$

**$^1\text{H}$  is 6400 times more sensible than  $^{13}\text{C}$**

Nuclei with larger  $\gamma$  will absorb/emit more energy, and will therefore be more sensitive.

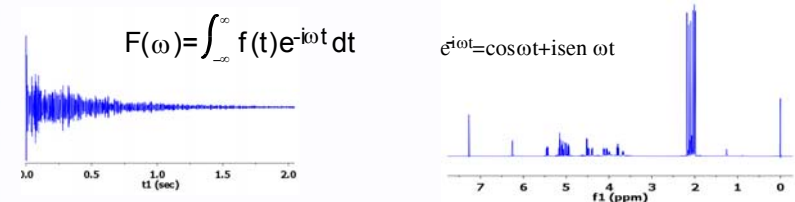
Unitat de RMN Centres Científics i Tecnològics

## Time Domain to Frequency Domain Fourier Transform


 Universitat de Barcelona




Time domain signals are converted into frequency domain signals using the **Fourier Transform**

$$F(\omega) = \int_{-\infty}^{\infty} f(t)e^{-i\omega t} dt$$

$$e^{i\omega t} = \cos\omega t + i\sin\omega t$$


*f(t)* corresponds to the time domain, and *F(ω)* corresponds to the frequency domain


*F(ω)* is a complex function that has a real (Re) and an imaginary part (Im)

Re	Absortion	A	B	C
Im	Dispersion			
		absortion signal	dispersion signal	absolute value signal $\sqrt{(Re)^2 + (Im)^2}$


line shape is Lorentzian (Fourier transform of a decaying exponential function)

Unitat de RMN Centres Científics i Tecnològics

## Frequency Domain

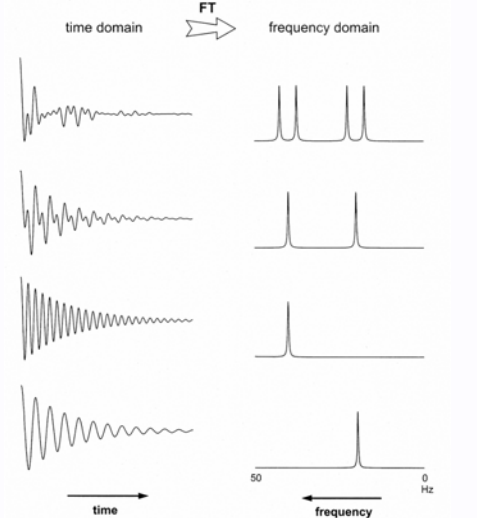
 Universitat de Barcelona

time domain  $\xrightarrow{FT}$  frequency domain



Unitat de RMN Centres Científics i Tecnològics

